


D₉

Anaesthesia



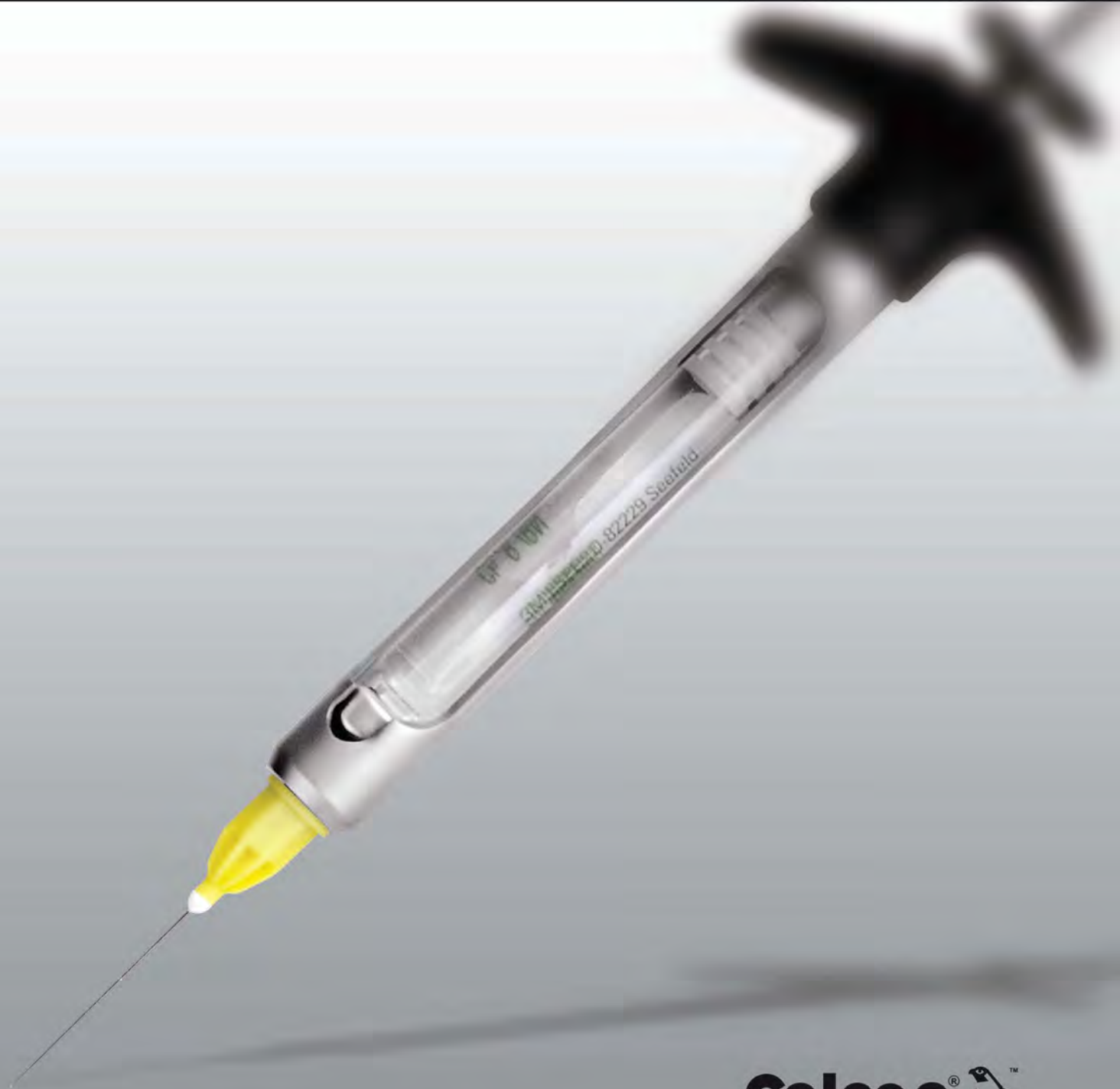

Stainless steel series
Made of 18-8 stainless steel
Longer life than Brass syringes
No chipping of chrome
Guaranteed quality
Sterilizable by any method

SS•2000

Stainless Steel Series

Intraligamental syringes



DA.010.180
Falcon-Ject Ultima
1,8 ml.

Injects 0.05 ml with every click
Supplied with angled and straight tubes



DA.040.180
Pistol-Ject
1,8 ml.

Injects 0.20 ml with every click
Lever action for high injection pressure

Dual aspiration system, self and manual aspirating syringes

During injection, aspiration can be achieved either manually by retracting the plunger backward or simply releasing pressure on plunger will produce self aspiration.



DA.626.180
Quatro
1,8 ml.



DA.726.180
Trio
1,8 ml.



DA.517.180
Type A
1,8 ml.



DA.557.180
Type C
1,8 ml.

Intraligamental syringes

Injects 0.06 ml with every click



DA.094.180
Swedish Type
1,8 ml.



DA.020.180
Falcon-Ject angled
1,8ml.



DA.021.180
Falcon-Ject straight
1,8ml.

Manual aspirating syringes



DA.801.180
Euro
Breech loading
1,8ml.



DA.806.180
Euro
Breech loading
1,8ml.



DA.635.180
T-Handle
1,8 ml.



DA.595.180
Euro
1,8 ml.



DA.705.180
3- Ring
1,8ml.



Made of Brass
 Chrome plated
 All Brass series syringes have stainless steel plungers
 Guaranteed quality
 Sterilizable by any method



Intraligamental syringes

Injects 0.06 ml with every click



DA.095.180
 Swedish Type
 1,8ml.

Injects 0.20 ml with every click
 Lever action for high injection pressure



DA.060.180
 Pistol-Ject
 1,8ml.



Plastic protection tube



DA.035.000
 Autoclavable

Disposable dental needles

Sterile, thread in compliance with European standard



- DA.1000.12 0,3 x 12 mm
- DA.1000.16 0,3 x 16 mm
- DA.1000.21 0,3 x 21 mm
- DA.1000.25 0,3 x 25 mm
- DA.1001.25 0,4 x 25 mm
- DA.1001.36 0,4 x 36 mm

Pack of 100 pieces

Manual aspirating syringes



DA.520.180
 Type A
 1,8ml.



DA.560.180
 Type C
 1,8ml.



DA.620.180
 Euro
 1,8ml.



DA.721.180
 3- Ring
 1,8ml.



DA.735.180
 2-Ring
 1,8ml.



DA.820.180
 Euro
 Breech loading
 1,8ml.




DA.821.180
 Euro
 Breech loading
 1,8ml.

P-Ject

Manual aspirating syringe

Special aspiration tip

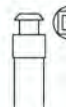


DA.561.180
 1,8 ml



Self aspirating syringes

Self aspiration: How it works

During use, the thumb pressure on the plunger or thumb disc stretches the diaphragm over a small metal projection inside of the syringe nozzle, increasing pressure inside of the cartridge. As the thumb pressure is released, pressure inside the cartridge lessens and aspiration occurs automatically. (Instruction for use are supplied with the product)

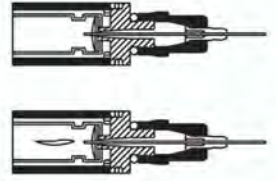

Disc system syringes

Dual self aspiration system:

Disc system syringe has dual aspiration feature to achieve self aspiration.

Aspiration before injection:

Insert syringe in tissue and press thumb disc to achieve aspiration.

Aspiration during injection:

Automatic aspiration occurs whenever pressure on the plunger is released.

SS•2000
Stainless Steel Series



DA.115.180
Disc system
1,8ml.

feather-lite



DA.120.180
Disc system
1,8ml.

Brass
SERIES



DA.145.180
Disc system
1,8ml.

Projection system syringes

Aspiration during injection:

Automatic aspiration occurs whenever pressure on the plunger is released.



SS•2000
Stainless Steel Series



DA.315.180
Projection system
1,8ml.

feather-lite



DA.320.180
Projection system
1,8ml.

Ultima

Needle re-capper

SS•2000
Stainless Steel Series

Falcon Needle Re-Cappers are the ideal way to prevent needle stick injuries when recapping needles. Place the upturned needle cap into Falcon Re-Capper and re-sheathing becomes a simple one hand procedure.



Needle re-cappers

SS•2000
Stainless Steel Series

Brass
SERIES



DA.980.000



DA.985.000



DA.990.000

